



made for you!



JOYCEone

JOYCEone

Your number 1 choice for a perfect fit in no time.

It has never been so quick and easy to fit a mask. The patient simply puts on the mask and the fitting is practically complete.

For all face shapes, design based on expertise with masks

Yesterday the patient had to try on mask after mask and go through a time-consuming fitting process. Today we offer a ventilation mask that fits nearly everyone. JOYCEone, the result of decades of experience and mask expertise, is made for you.

Hand-in-hand: comfort and effective therapy

Thanks to its optimized shape, the single-layer mask cushion finds its ideal position, regardless of the size and shape of the patient's nose. The forehead support works just as cleverly. The support's springs automatically find and maintain the right distance to the forehead to ensure effective therapy and patient comfort.

With attention to detail

Your patients will fall in love with JOYCEone at first sight! This mask boasts more than a perfect fit and wearing comfort. With the unmistakable color coding of headgear and clips, your patient can assemble JOYCEone correctly and quickly right from the start. What's more, cleaning the mask is more convenient and easier than ever.


Modern Efficiency

Short and sweet: the intelligent JOYCEone helps you to save time without compromising on patient care. JOYCEone – made for you!

JOYCEone Advantages

- Economical: one size fits nearly all patients
- Time-saving perfect fit: automatic adjustment of mask cushion 1 and forehead support 3
- Unmistakable assembly: color coding on headgear and clips 4
- Convenient: effortless cleaning
- Pleasant: quiet exhalation system 7
- Stable: very durable materials

Noticeably greater wearing comfort and assured cleanliness: Our coating WEINMANN SilkTec of liquid silicone rubber sets new standards with its significantly smoother and more hygienic surface than conventional coatings.

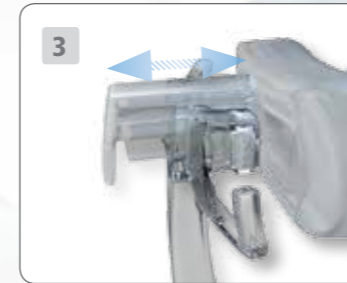






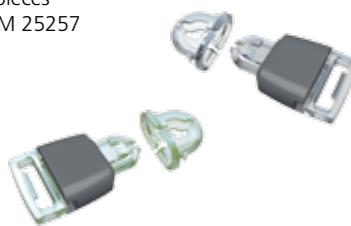



- 1 Mask cushion
- 2 Forehead cushion
- 3 Forehead support
- 4 Headgear clips
- 5 Hook-and-loop tab
- 6 Ball-and-socket joint
- 7 Exhalation system
- 8 Pressure measurement port




JOYCEone

Product name	WM no.
JOYCEone	WM 25260


Replacement parts

Headgear incl. headgear clips WM 25255 	Forehead cushion silicone, WEINMANN SilkTec coating WM 26801 	Connection for pressure measurement port WM 26849 
	Forehead support incl. silicone forehead cushion, WEINMANN SilkTec coating WM 15907 	Rotating sleeve WM 26846 
Headgear clips 2 pieces WM 25257 	Mask cushion silicone, WEINMANN SilkTec coating WM 25253 	Ball-and-socket set, vented retaining ring, elbow, rotating sleeve, plug WM 15908 
		Carrying bag WM 26295 

Our complete selection of therapy solutions, accessories and mask systems can be found at: weinmann.de

Technical Data JOYCEone		CE 0197 	
Product class as per Directive 93 / 42 / EWG	II a	Temperature range	<ul style="list-style-type: none"> Operation: +5 °C to +40 °C Storage: -20 °C to +70 °C
Dimensions (W x H x D)	ca. 96 x 115 x 83 mm	Flow resistance	<ul style="list-style-type: none"> at 50 l/min: 0.05 hPa at 100 l/min: 0.28 hPa
Weight	87 g	Quoted two-figure noise emission value as per ISO 4871	<ul style="list-style-type: none"> Sound pressure level: 15 dB(A) Sound power level: 23 dB(A) Uncertainty factor: 3 dB(A)
Therapy pressure range	4 hPa bis 30 hPa	Service life	up to 12 months ¹
Tube connection Tapered connection as per EN ISO 5356-1	Ø 22 mm (male)		
Pressure measurement port	Ø 4 mm		

¹ Materials used in the manufacture of masks age when exposed to aggressive cleaning products. Therefore, it may sometimes be necessary to replace mask parts earlier.

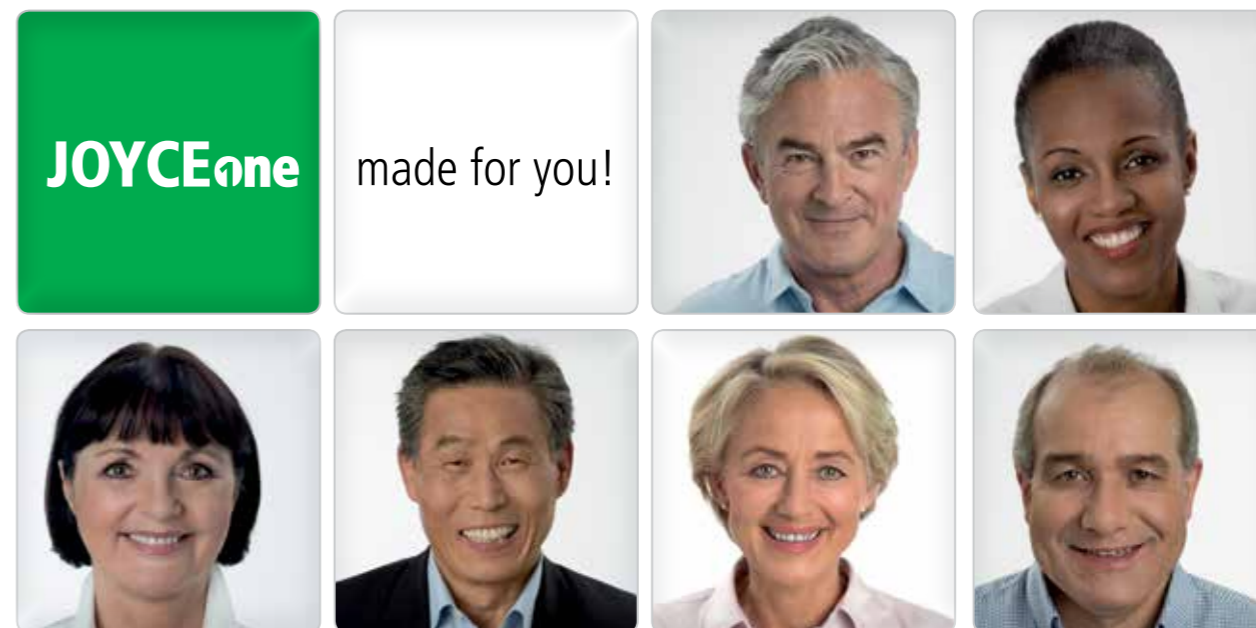

JOYCEone

83699-EN-08/2013 © Copyright Weinnmann, Hamburg. Duplication of any kind only with the express permission of Weinnmann. All rights to design and specification modifications reserved.